

sagemax.com

sagemax®



NexxZr. Glaze Spray

Effective glazing.  
Efficient processing.

Effective glazing.  
Efficient processing.

# NexxZr® Glaze Spray



**NexxZr® Glaze Spray is a spray glaze made of extremely high-strength lithium silicate glass-ceramic for the finalization of monolithic zirconium oxide restorations.**

The optimal size of the glass particles as well as the finely adjusted spray nozzle enable a light, even and very thin application of the material. Already one thin layer provides for a high adhesive bond to zirconium oxide and results in a smooth, homogeneous surface - preserving even the most intricate surface details. The antagonist-friendly wear characteristics of the lithium silicate glass-ceramic provide a long-lasting protection for high-stress areas (e.g. contact points).

## Technical data



**Type**  
Lithium silicate  
glass-ceramic



**Flexural strength**  
> 140 MPa



**Light transmission**  
> 90%, clear



**Bond strength**  
> 35 MPa



**Abrasion**  
Not measurable

## Highlights

---

1

### High-strength lithium silicate glass-ceramic

High adhesive bond with  $ZrO_2$  and smooth, homogeneous surfaces



Extremely thin material application approx.  $10\ \mu\text{m}$  and high adhesive bond (conventional glaze pastes approx.  $40\text{--}50\ \mu\text{m}$ )


Smooth, homogeneous surfaces without loss of surface structure

---

2

### Optimal glass particle size and finely adjusted spray nozzle

Light, even and very thin material application in seconds

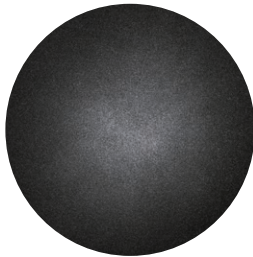

Extremely fine spray



Light, uniform and very thin material application in seconds

---


3

### Antagonist-friendly wear resistance

Long-lasting protection for high-stress areas, e.g. contact points



High abrasion resistance for long-lasting protection

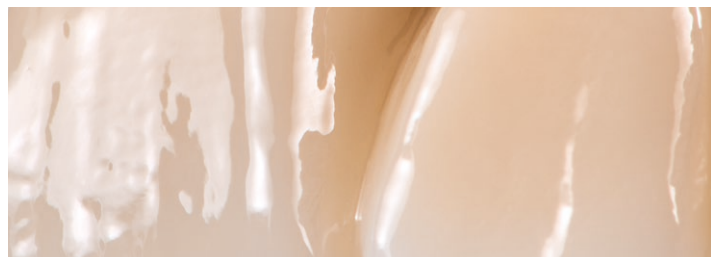

Smooth, homogeneous surfaces

## Time & Consumption



\* Internal test

## Delivery forms

**Spray bottle**  
75 ml

## Other top products



**NexxZr+ Multi**  
A multilayer zirconium oxide with natural esthetic



**NexxZr T Multi**  
A multilayer zirconium oxide with multifunctional esthetics and a wide variety of indications



**NexxZr+**  
A zirconium oxide with high translucency for esthetic anterior and posterior restorations



**NexxZr T**  
A zirconium oxide with medium translucency and high flexural strength for a variety of indications and processing options

**Sales & Technical Support EMEA**  
EC Representative AB Ardent  
Generatorgatan 8, 195 60 Arlandastad, Sweden  
T: +46 8 594 412 57, E: info@sagemax.com  
sagemax.com